


Build your Dispersing
system by selecting a
motor plus one Tool

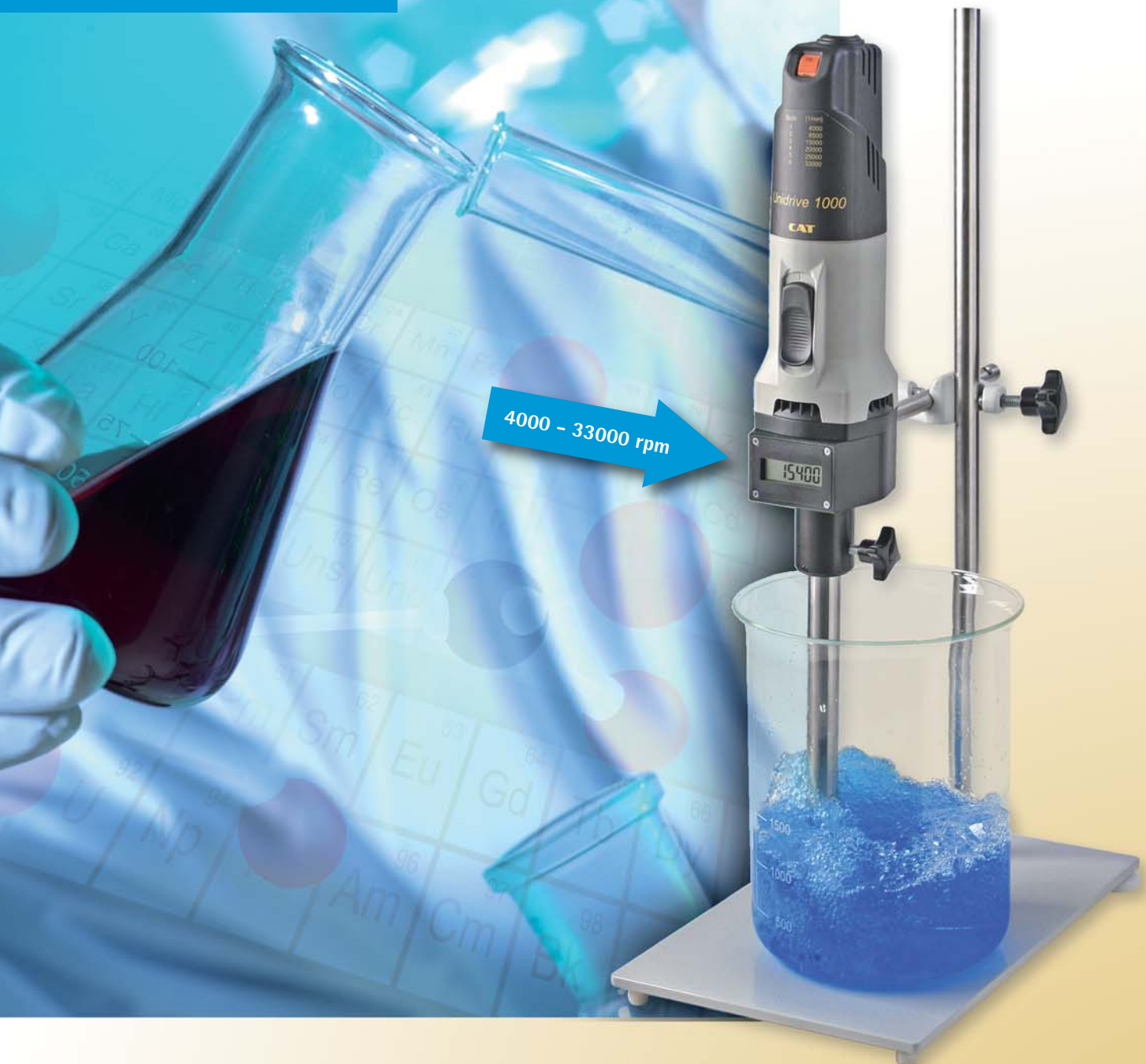

4000 - 33000 rpm

From marketing knowledge to real innovation

Homogenisers

X1740



5.000 - 20.000 ml
2.500 - 23.000 rpm

T40 - 20.000 ml (CAT No.: 60415)

X1740

high-power, heavy-duty drive motor with electronic speed control, 1700 watts for 40 mm tool 60408-00

0,2ml - 5.000ml
4.000 - 33.000 rpm

Unidrive

powerful drive motor, built-in speed control for all dispersing tools below up to 30 mm 60402-50

T20 250 ml - 2.000 ml (CAT No.: 60423)

T30 500 ml - 5.000 ml (CAT No.: 60425)

Unidrive D

same as above but with digital display of rpm 60406-50

0,2ml - 2.000ml
45.000 rpm

T20 250 ml - 2.000 ml (CAT No.: 60423)

X360

drive motor, 350 watts, high speed motor, speed control separately available 60401-00

X120

hand-held, high speed drive motor for small volumes electronically stabilized speed control for tools up to 17 mm 60404-00

T20/G20	T30/G30	T40/G40

T = Shaft with ETFE bearing. This shaft type is used for standard application.

G = Shaft with slide seal ring. This shaft type is excellent for applications with aggressive respectively abrasive media and when pressure is applied.



Product Overview



Find out more about us,
our business fields and products under